Bedienungsanleitung MA2602 MA2603

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Nennzeit: 15/30min

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Instruction Manual MA2602 MA2603

- Please keep this product manual properly for easy reference at any time and carefully read it before using the product.

Directory

1、Catalogue	1
2、 Safety precautions	2
(1) Usage environment	2
(2) Those who are not suitable for using this device	2
3、 Preparation before use	3
(1) Seat cushion installation	3
(2) Sitting posture adjustment	3
(3) Method of Moving the Massage Chair	4
(4) Placement of massage chairs	4
(5) Starting power supply	5
4、 Main construction names and functions	6
5、 Introduction to Control Panel Functions	7
6、 Massage Program Operating Instructions	8
7、 Daily maintenance	9
8、 Technical parameters	10

Reminder:

Our company reserves the right to make improvements to the product and instructions without prior notice. The appearance of the product shall be based on the actual product. If there are differences between the standard configuration description and the actual product due to product upgrades or other reasons, please refer to the actual product.

Safety precautions

(1) Usage environment

- (1) After use or before cleaning the product, please unplug the power plug from the socket.
- (2) Do not use or store this device in particularly humid environments such as bathrooms, saunas, hot springs, etc. to avoid mold, electric shock, fire, or mechanical failure
- (3) Do not use in environments with room temperature above 40 degrees Celsius
- (4) Do not use this device outdoors.
- (5) Do not use this machine in an environment with too much dust, oil stains, or oxygen deficiency. (6) Do not place this machine in a hot air stove, stove, or direct sunlight, as this can cause overheating of the parts inside the massage chair and cause fire.
- (7) Do not use this machine in an uneven area. Please place the support in a level area such as the floor for use.
- (8) To prevent electric shock, please confirm whether the socket is grounded. Please consult the sales store when installing the grounding wire for the socket
- (9). Do not use it while smoking.

(2) Those who are not suitable for using this device

The following personnel should consult a doctor before using (which may cause physical discomfort or even accidents, resulting in danger).



- (1) Patients with malignant tumors or malignant pus masses.
- (2) Patients with severe heart disease.
- (3) Patients with severe skin diseases.
- (4) Patients with severe bone hyperplasia.
- (5) Huai Zhe.
- (6) Under the guidance of a doctor for recuperation, especially for those who feel physically abnormal
- (7) Individuals with spinal abnormalities or back curvature caused by accidents or illnesses.
- (8) Those who need quietness.
- (9) Joint deformers.
- (10) Individuals with medical devices such as pacemakers implanted in their bodies.
- (11) Individuals with acute symptoms of visceral diseases such as gastritis, enteritis, dementia, hepatitis, etc.
- (12) Individuals with severe blood circulation disorders in the legs.
- (13) Individuals with sensory impairments must pay attention when using this product as its surface heats up and is not sensitive to heat.
- (14) Individuals with high fever symptoms.
- (15) This product is not recommended for use by individuals (including children) who have physical disabilities, sensory or mental impairments, or lack relevant experience and knowledge, unless there is supervision and guidance from someone responsible for their safety when using this product.
- (16) Children should be taken care of to ensure that they do not play with this product.

Do not allow individuals under the age of 18 or those who cannot control their behavior (such as mentally ill individuals, mentally disabled individuals, paralyzed patients, etc.) to operate the massage chair alone or come into contact with various parts of the chair. This situation may lead to electric shock, injury to the user, or component failure.

Preparation before use


(1) Seat cushion installation

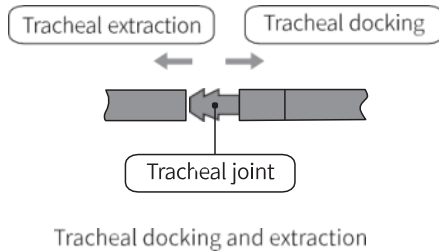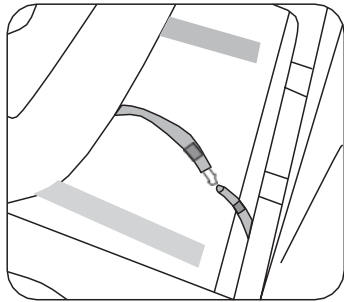(2) Seat cushion adjustment

• Press the relevant buttons such as "backrest up" and "backrest down" to adjust to the most comfortable angle for the body

※When pressing and holding the chair back lift button, the chair back (feet) will continue to rise or fall. When tilted to a suitable position, simply release the button to stop.

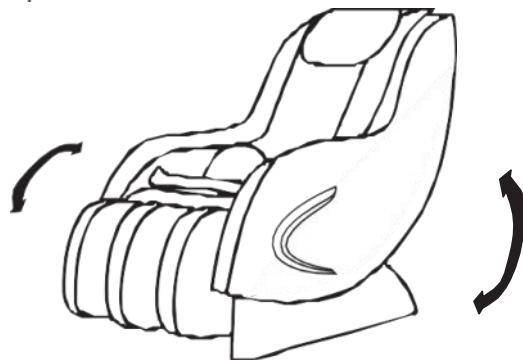Requirement: When adjusting the backrest and leg pads, ensure that there is sufficient space left in the front and back of the massage chair.

(3) Method of Moving the Massage Chair

When moving, lift the handrails on both sides and move.

Attention!


Do not move with people.

Do not use small pulleys to move on easily damaged floors.

(4) Placement of massage chairs

Please leave enough space around the massage chair when preparing for the massage


- (1) The place where the massage chair is placed should have sufficient space for the chair back to recline and the foot rest to rise.
- (2) Ensure a space of at least 20 centimeters at the back and at least 20 centimeters in front of the massage chair.
- (3) Ensure that there are no obstacles in front and behind the massage chair.
- (4) To prevent damage to the ground, please place a carpet under the massage chair to avoid damaging the ground.

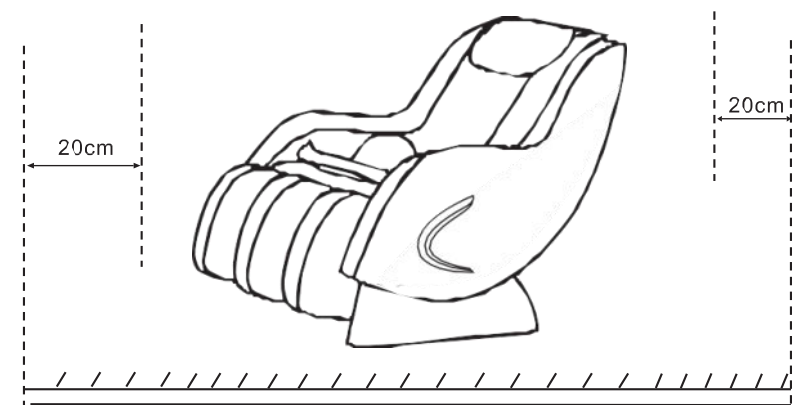(5) Starting power supply

(1) Confirm the surrounding environment




Confirm that there are no obstacles around the massage chair

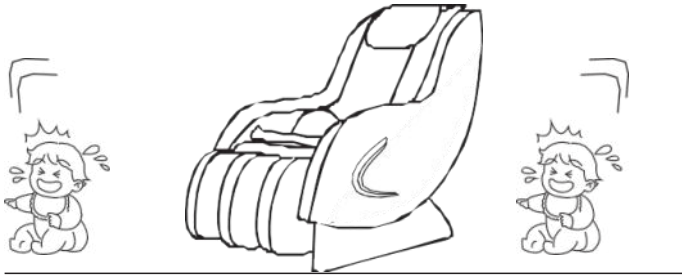
(2) Insert the plug into the power socket and turn on the power switch

1. Please check and confirm if the power cord is pressed under the massage chair
2. Is the ground wire properly connected
3. Before inserting the plug, ensure that all operating buttons (including the main power switch) are turned off

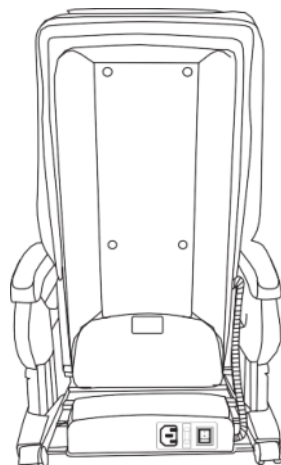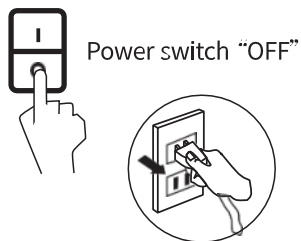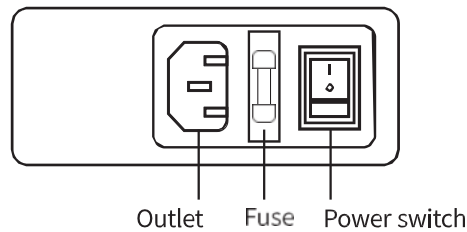Main construction names and functions

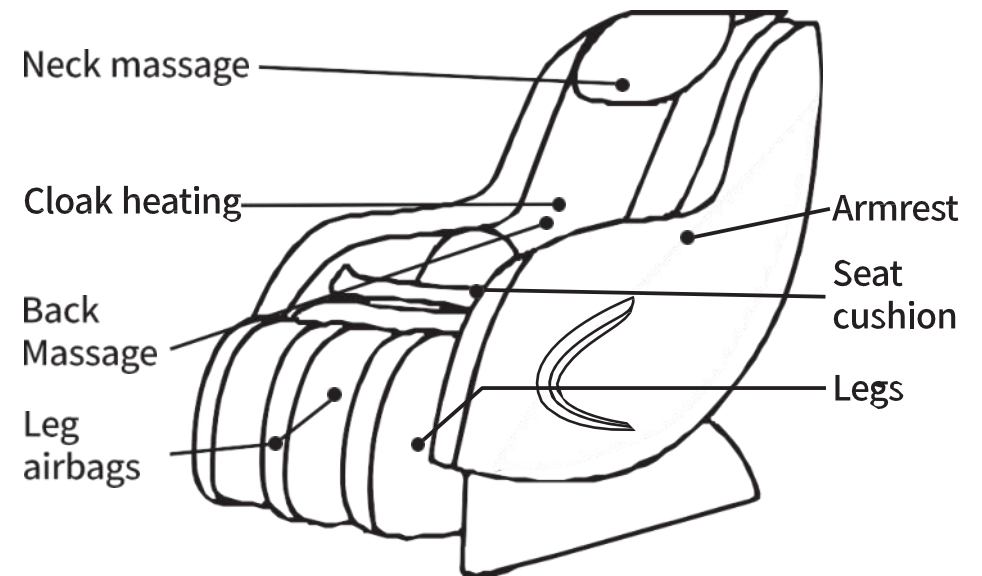Instructions for operating buttons

Start button: defaults to turning on the machine and entering 15 minute automatic mode. Press again to enter 30 minute mode, and press again to turn off the massage chair


Mode: There are multiple automatic modes to choose from

Intensity: Adjust the intensity of the back massage, with three adjustable gears: high, medium, and low

Heating: Turn heating function on/off

Rise: Long press "Rise" to raise the backrest

Drop: Long press "Drop" to lower the backrest

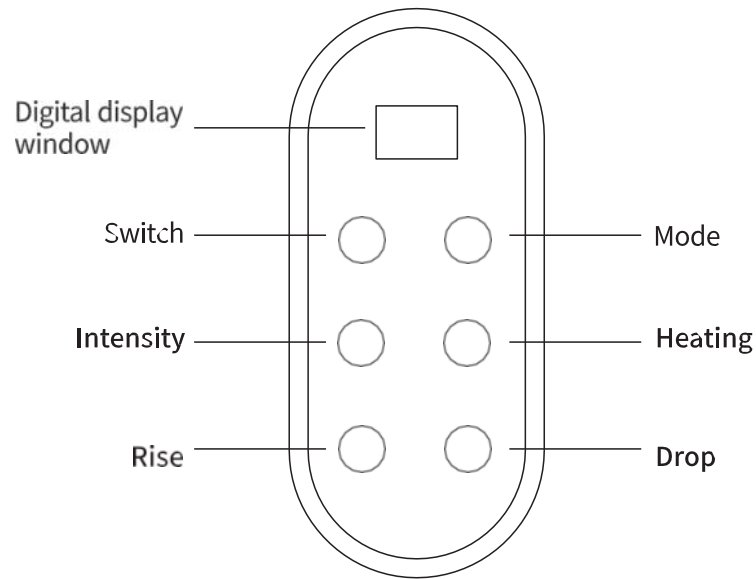Instructions for operating buttons

On key: By default, the machine will enter 15 minute automatic mode when turned on, and then press again to enter 30 minute mode. Press again to turn off the massage chair;

Mode : There are multiple automatic modes to choose from;

Intensity: Adjust the intensity of the back massage, with three adjustable gears: high, medium, and low.

Heating: Turn on/off the heating function

Rise: Long press "Rise" to raise the backrest

Drop: Long press "Drop" Backrest down

Product parameter details

The entire machine adopts blow molded handrails; The back adopts a suspended massage chair movement that automatically adjusts the angle according to different body types. One piece installation free design, ready to use after unpacking; Arbitrarily adjusting the angles of the calves and backrest to achieve zero gravity massage of the spacecraft; Metal movement, silent design, using PU massage chair specialized leather material, and exquisite sewing technology; Additional features include a built-in Bluetooth speaker to connect to your phone and play music.

Product features

Massage areas: Head, Neck, Back, Waist, Legs

Airbag part: Leg

Massage techniques: Kneading, Massage, Squeezing, Hot compress

Massage intensity: You can adjust the massage intensity according to your own situation

Legs: Airbag massage legs

Daily maintenance

- (1) The massage chair should be kept clean and dust-free at all times, and placed in a non damp place
- (2) If not used for a long time, please cover it with a cloth.
- (3) Do not place the massage chair in direct sunlight or at high temperatures, as this can cause the leather of the massage chair to fade and accelerate aging.
- (4) To clean the massage chair, please turn off the main power switch and unplug the power plug
- (5) Do not use cleaning agents containing organic solutions, such as benzene, to clean massage chairs, as this may cause leather discoloration, aging, cracking, etc
- (6) If the power cord is damaged, it must be replaced with a dedicated cord or a dedicated component purchased from its manufacturer or maintenance department.

Technical parameters

Product model: MA2602 MA2603 Size: 140cm

Rated voltage: 110V-240V

Net weight/gross weight: 41kg

Rated power: 70W

Standing size: 119 * 72 * 105cm

Rated time: 15/30min

Applicable height: 155-178cm

Safety level: Class I